

serena[®]
pillow lift


Operating Instructions

 mountway[®]
bringing independence to life

Please read these instructions thoroughly before using your Serena Pillow Lift.

1. Contents

- Serena Pillow Lift
- Serena mattress and cover
- Control box
- Hand control

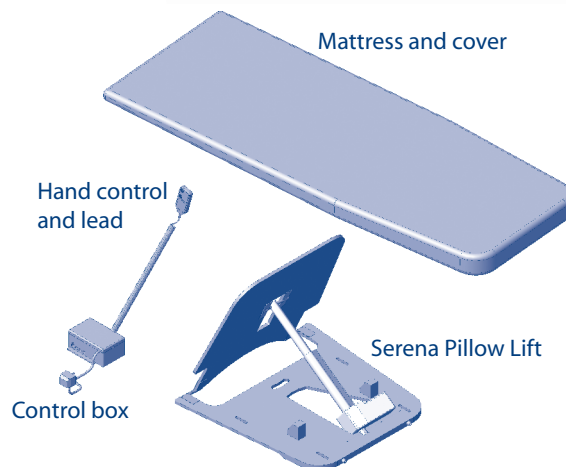

Please inform our customer services department on 01495 723300 within seven days of receipt if any contents are damaged or missing.

2. Using the Serena

- Place Serena on your bed, clear of the headboard
- Plug the lead from the hand control into the control box
- Plug the lead from Serena's actuator into the control box
- Place the control box on the floor in a convenient and safe position and plug into a mains power point
- To raise Serena, press the UP button on the hand control. Release this button when Serena reaches the desired position
- To lower Serena, press the DOWN button

3. Cleaning

- **Frame** - wipe clean with a non-abrasive household cleaner and dry thoroughly
- **Mattress cover** - wipe clean with a non-abrasive household cleaner or hand/machine wash to 95°C. Dry by hanging out, spinning or tumble drying at temperatures up to 130°C
- **Mattress** - this cannot be washed and should be replaced if soiled. Call our customer services department on +44 (0)1495 723300 for information

4. Precautions

- Check that the operating parts of Serena are free of obstructions before and during use
- Do not use Serena's components with any other products
- Serena contains no user serviceable parts. In the unlikely event of a fault, please contact our customer services department on +44 (0)1495 723300
- If Serena has been moved from the bed e.g. for changing sheets, on replacing make sure all electrical connections are securely fitted to the control box

Mountway Limited

Tafarnaubach Ind. Est. Tredegar NP22 3AA UK

tel +44 (0)1495 723300 fax +44 (0)1495 723360

e-mail mail@mountway.co.uk

www.mountway.com

